

μDrop GPS® Dish

Embryo Culture Dish for Human IVF

Safe

Effective

Easy-to-use



Jacques Cohen, PhD

"The Micro GPS dish evolved from the input of dozens of scientists and engineers. It is the most versatile and ingenious culture dish ever developed. It combines so many new features.

It is the strongest and safest Petri dish for culturing embryos. The bottom underneath each GPS location is so thin, you can see the embryo smile."




μDrop GPS®

New Design Features

- Precise 20 μl micro-wells with GPS® feature for rapid location and visualization
- Enhanced optics
- Better orientation and identification
- EmbryoAddress™ helps to quickly identify and track embryos
- 32 mm inner 'oil ring' for VOC protection
- Uses up to 75% less oil than a conventional 60 mm dish
- Fill the outer wells with oil for better temperature stability out of the incubator
- Works well with both oil-overlay and oil-underlay methods
- Designed to ensure safe embryo culturing
- For use with both standard and mini incubators

The working surface and well bottoms

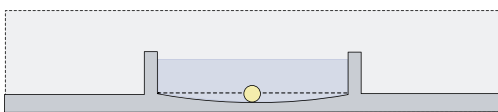

● = Egg or Embryo

New Design Features

The working surfaces (bottom) of all wells of the μDrop GPS® dish is sloped to create the GPS® location for all the specimens within these wells. This makes locating them quicker, easier to handle and observe, and the ability to return specimens to their controlled environment sooner.

General Features of GPS® Dishware

- Improved safety (no droplet collapsing or mixing)
- Better pH control
- GPS® micro-wells protect against sudden movements
- Raised lid promotes gas exchange, prevents contamination, and stops oil seals from forming
- Non-toxic medical grade, non-pyrogenic polystyrene
- Better ergonomics than a 35 mm dish
- CE and ISO Registered, FDA 510(k) Cleared
- 1-cell MEA and < 0.03 LAL Tested

New Breathable Packaging for all GPS® Dishware

- Reduces off-gassing time
- Less VOCs introduced into the laboratory
- Tested and validated to maintain sterility for the entire 5-year shelf life of the dish

IVFonline.com®

 **LifeGlobal®**
GROUP

T: 1-800-720-6375

F: 1-519-826-6947

Intl.: 001-519-826-5800

sales@LifeGlobal.com

www.LifeGlobalGroup.com

©2015 LifeGlobal® Group, LLC

Dish_uDrop-2015-12-01-1155